
USER'S GUIDE AGATA ANALYSIS

SOFT CUBIX

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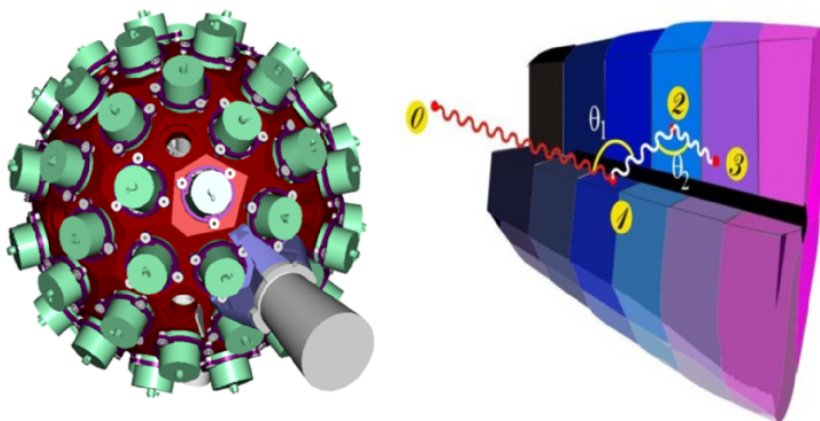```

** Print this help :  $p + h$ 
** Histogram selection ** ==> Click on a histogram make it the active histogram (the default one is the first
histogram in the list of primitives of the active pad)
** Histogram operations **
==>  $b + e$  : background evaluation ==>  $b + s$  : background subtraction
** background operation **
==>  $i + b$  : Starting interactive background adjustment utility
If the background adjustment utility is on:
==>  $s + i$  : active the iterator variations by mouse wheel
==>  $s + f$  : active the filter order variations by mouse wheel
==>  $s + s$  : active the Smoothing windowd variations by mouse wheel
==>  $s + c$  : active/unactive compton
==>  $s + d$  : change direction (BackDecreasingWindow/BackIncreasingWindow) Press Escape key to exit the utility
** Peaks operations **
==>  $h + m$  : Hide markers
==>  $s + m$  : Show markers
==>  $a + s$  : Adding a gamma source from the nuclear database
==>  $f + a$  : Open fit all panel menu ==>  $a + p$  : Add a peak at mouse position
==>  $A + P$  : Add a peak + background at mouse position
==> CTRL (  $A + P$  ) : Add a point to the graphical polyline
==>  $s + l$  : Save the peaks of the current pad the database
==>  $l + l$  : Load a peak list from the data base ==>  $c + a$  : Clear all (peaks + graphical polyline)
==>  $o + m$  : Open PeakCreator menu
==>  $c + 0$  : Switch Color to default #0
==>  $c + 1$  : Switch Color to default #0
==>  $c + 2$  : Switch Color to default #0
==>  $c + 3$  : Switch Color to default #0

```

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```
TString Raw_Directory="/agataDAS1/AGATA/Experiments/e680/Analysis/Users/Gui/Level_Scheme/Import_ENSDF/";
TString Tree_Directory="/agataDAS1/AGATA/Experiments/e680/Analysis/Data/AnalysedTrees/E680_PerIsotope_PerZ_PerM_NoSelection.root";
const char* Exp_Directory = "/agataDAS1/AGATA/Experiments/e680/Analysis/Users/Gui/Level_Scheme/Import_ENSDF/EXP_DATA/";
const char * E_AGATA_DIRECTORY="/agataDAS1/AGATA/Experiments/e680/Analysis/Users/Dudouet/Calibrations/BuildTrees/Trees/Runs.186_188.root";
```

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```
*****
*****HELP AND PROGRAM Last UPDATE : 19/06/2016 15h16*****
*****
HELP : Keys used to interact with nucleide chart:
h->Show Help
CTRL+ h->Show Special Help for Gated Spectra
r->Reset coordinates of the chart
z->Change neutron by proton
s->Change proton by neutron
d->Increase by one neutron
q->Increase by one proton
Ctrl+a -> Draw All Nuclei
Ctrl+r -> Reset Drawing

-----
Ctrl+d -> Draw Spectrum of one Nuclei ex: 100Zr. Option : L(Draw LS), V(Verbose Mode)
OPTION2(a) : Z+1;Z-1;Z+- -> Draw Energies of the Z+(-)(+and-) Nuclei
OPTION2(b) : A+1;A-1;A+- -> Draw Energies of the A+(-)(+and-) Nuclei
OPTION2(c) : A++;A-;A++- -> Draw Energies of the A+2(-2)(+2and-2) Nuclei
OPTION2(d) : A+++;A-;A+++- -> Draw Energies of the A+3(-3)(+3and-3) Nuclei
OPTION2(e) : A+4;A-4;A++++- -> Draw Energies of the A+4(-4)(+4and-4) Nuclei
All Options (a),(b),(c) can be called together separated with comma
Ex: if one want display all contaminants, he should write :
Z+- A+- A++- A+++- A++++-
Accepted forms : 'Z+1','z+1','Z+','z+', 'Z-1','z-1','Z-','z-', 'Z+-','Z-+', 'z-+'
Accepted forms : 'A+1','a+1','A+','a+', 'A-1','a-1','A-','a-', 'A+-','a+-', 'A-+', 'a-+'
Accepted forms : 'A+2','a+2','A++','a++', 'A-2','a-2','A-','a-', 'A++-','a++-', 'A-+-','a-+-'
Accepted forms : 'A+3','a+3','A+++','a+++', 'A-3','a-3','A-','a-', 'A+++-','a+++-', 'A-+-+', 'a-+-+'
Accepted forms : 'A+4','a+4','A++++','a++++', 'A-4','a-4','A-','a-', 'A++++-','a++++-', 'A-++++', 'a-++++'
-++++'
StrengthMin : Display only intensities above this limit
SAME : Display new nuclei on the same Pad, one can also displays energies of new contaminants
Quality_Factor : Trust Factor of the particle identification
Accepted forms : 'X/Y' Where X is Z_ID and Y M_ID in case of gated/raw spectra per isotope
Accepted forms : 'X' or 'X/Y' Where X is M_ID in case of gated/raw per Mass and Z_ID/M_ID in case of E
versus M/Z Matrix. Y is not taken into account.
GATE : Set of gates ( see the special Help )
GATE_W : Gate Width
Max_Entries : Limit the entries readed

-----
Ctrl+i -> Draw Spectra of an isotopic set between Amin and Amax
Ctrl+w -> Draw Distribution of possible pollutants weighted by probabilities
Ctrl+e -> Draw Level Scheme determined in this work.
All the GLS keyboard shortcuts are available to modify the Level scheme
You can save the LS in the current directory, and it will be reused to draw arrow of possible contaminants
Ctrl+k -> Connect the canvas selected in order to use special BackGround subtraction
Ctrl+b -> Draw Theoretical raw intensities determined in this work by Loading a LvlScheme
Energy_Ref -> Energy used as a reference in order to display other intensities
Height -> Height of the Energy reference peak
Ctrl+n -> Draw Theoretical gated intensities determined in this work by Loading a LvlScheme !!!! CTRL+n
DONT'T WORK FOR THE MOMENT 12/02/2016
Ctrl+x -> Exit
```

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CTRL+ h->Show Special Help for Gated Spectra

HELP ON OPTIONS TO PLOT GATED SPECTRA

KEYWORDS :

Default :N 0 -> No Gates

SM-> Put simple Gate on Tree conditioned by Mass

SI-> Put simple Gate on Tree conditioned by Mass and Z

DM-> Put Double Gate on Tree conditioned by Mass

T-> Put Triple Gate ...

SA-> Put Simple Gate AGATA ONLY

SYNTAX :

SM 125+256+265

One gate with several 'or' conditions

DM 125.3+659.3/254.5+256.9

Double Gate with 'or' connections

HELP ON OPTIONS TO PLOT ISOBAR OR ISTONES MATRIX

E_Z-> Plot E versus Z Matrix with current M between Current Z-3 and Current Z+3

E_M-> Plot E versus M Matrix with current Z between Current M-3 and Current M+3

HELP ON OPTIONS TO PLOT RAW SPECTRA

PIDM-> Plot PID versus M for the current Mass

B_M-> Plot Raw Spectra with the current mass

Options are :

B_M->Raw Spectra with all the available PIDs

B_M X+/-x-> Raw Spectra for the current mass and PID X+/-x

HELP ON OPTIONS TO Subtract And Mutliply/Divide by another Spectra

S_S-> KeyWord To substract Current Spectra by Another

Options are :

S_S X N 0->Add to the Current Spectra another Spectra Corrected by factor X

S_S -X SI 369+542->Subtract to the Current Spectra the Gated Spectra of the current Nuclei Gated on 369+542

S_S -X DM 369+542/586->Subtract to the Current Spectra the Double_Gated Spectra of the current Nuclei

HELP ON OPTIONS TO SEARCH FOR CONTAMINANTS

S_P-> KeyWord To search for contaminants

Options are :

S_P W+X/Y+Z -> W and X are Z_ID and M_ID for the fisrt Spectrum; Y and Z are Z_ID and M_ID for the second Spectrum

Accepted forms : 'S_P 0.5+0.5/0.321+0.321'

Accepted forms : 'S_P 0.321+0.321' -> The second parameters are taken from teh default values for M_ID and Z_ID

Default Values: 'Z_ID=M_ID=0.5'

HELP ON OPTIONS TO USE AUTOMATIC BACKGROUND SUBTRACTION

ES-> KeyWord to use background subtraction

USE :

SM X+Y ES

ES is put at the end of the gate expression

Warnning! Work Only With the SM Mode